



Standard Test Method for Anions in Water by Chemically Suppressed Ion Chromatography¹

This standard is issued under the fixed designation D 4327; the number immediately following the designation indicates the year of original adoption or, in the case of revision, the year of last revision. A number in parentheses indicates the year of last reapproval. A superscript epsilon (ϵ) indicates an editorial change since the last revision or reapproval.

1. Scope

1.1 This test method² covers the sequential determination of fluoride, chloride, nitrite, *ortho*-phosphate, bromide, nitrate, and sulfate ions in water by chemically suppressed ion chromatography.

NOTE 1—Order of elution is dependent upon the column used; see Fig. 1.

1.2 This test method is applicable to drinking and waste waters. The ranges tested for this test method for each anion were as follows (measured in mg/L):

Fluoride	0.26 to 8.49
Chloride	0.78 to 26.0
Nitrite-N	0.36 to 12.0
Bromide	0.63 to 21.0
Nitrate-N	0.42 to 14.0
<i>o</i> -Phosphate	0.69 to 23.1
Sulfate	2.85 to 95.0

1.3 It is the user's responsibility to ensure the validity of this test method for other matrices.

1.4 Concentrations as low as 0.01 mg/L were determined depending upon the anions to be quantitated, in single laboratory work. Utilizing a 50- μ L sample volume loop and a sensitivity of 3 μ S/cm full scale, the approximate detection limits shown in Table 1 can be achieved. If lower detection levels are required, the sensitivity may be improved by using a lower scale setting (<3 μ S/cm) or a larger sample injection loop (>100 μ L). The analyst must assure optimum instrument performance to maintain a stable baseline at more sensitive conductivity full-scale settings.

1.5 The upper limit of this test method is dependent upon

total anion concentration and may be determined experimentally as described in Annex A1. These limits may be extended by appropriate dilution or by use of a smaller injection volume.

1.6 Using alternate analytical column and eluents may permit additional anions such as formate or citrate to be determined. This is not the subject of this test method.

1.7 *This standard does not purport to address all of the safety concerns, if any, associated with its use. It is the responsibility of the user of this standard to establish appropriate safety and health practices and determine the applicability of regulatory limitations prior to use.*

2. Referenced Documents

2.1 ASTM Standards:

- D 1066 Practice for Sampling Steam³
- D 1129 Terminology Relating to Water³
- D 1193 Specification for Reagent Water³
- D 2777 Practice for Determination of Precision and Bias of Applicable Methods of Committee D-19 on Water³
- D 3370 Practices for Sampling Water from Closed Conduits³

3. Terminology

3.1 *Definitions*—For definitions of terms used in this test method, refer to Terminology D 1129.

3.2 Definitions of Terms Specific to This Standard:

3.2.1 *analytical column set*—a combination of one or more guard columns followed by one or more analytical columns used to separate the ions of interest. It should be remembered that all of the columns in series contribute to the overall capacity of the analytical column set.

3.2.2 *chemical suppressor device*—a device that is placed between the analytical columns and the detector. Its purpose is to inhibit detector response to the ionic constituents in the eluent, so as to lower the detector background and at the same time enhance detector response to the ions of interest.

3.2.3 *eluent*—the ionic mobile phase used to transport the sample through the system.

3.2.4 *guard column*—a column used before the analytical column to protect it from contaminants, such as particulate matter or irreversibly retained materials.

3.2.5 *ion chromatography*—a form of liquid chromatography in which ionic constituents are separated by ion exchange followed by a suitable detection means.

3.2.6 *resolution*—the ability of an analytical column to separate constituents under specific test conditions.

3.2.7 *analytical column*—the ion exchange column used

¹ This test method is under the jurisdiction of ASTM Committee D-19 on Water and is the direct responsibility of Subcommittee D19.05 on Inorganic Constituents in Water.

Current edition approved Aug. 10, 1997. Published January 1998. Originally published as D 4327 – 84. Last previous edition D 4327 – 91.

² The following references may be consulted for additional information:

Small, H., Stevens, T. S., and Bauman, W. C., "Novel Ion Exchange Chromatographic Method Using Conductometric Detection," *Analytical Chemistry*, Vol 47, 1975, p. 1801.

Stevens, T. S., Turkelson, V. T., and Alve, W. R., "Determination of Anions in Boiler Blow Down Water with Ion Chromatography," *Analytical Chemistry*, Vol 49, 1977, p. 1176.

Sawicki, E., Mulik, J. D., and Witgenstein, E., Editors, *Ion Chromatographic Analysis of Environmental Pollutants*, Ann Arbor Science Publishers, Ann Arbor, MI, 1978.

Mulik, J. D., and Sawicki, E., Editors, *Ion Chromatographic Analysis of Environmental Pollutants*, Vol/No. 2, Ann Arbor Science Publishers, Ann Arbor, MI, 1979.

Weiss, J., *Handbook of Ion Chromatography*, Dionex Corp., Sunnyvale, CA, 1986.

Waters Innovative Methods for Anion Analysis, Waters Chromatography Division of Millipore, Method A-107 and A-116, 1990.

Haddad, P. R., and Jackson, P. E., *Ion Chromatography: Principles and Applications*, Elsevier Scientific Publishing Co., 1990.

³ *Annual Book of ASTM Standards*, Vol 11.01.

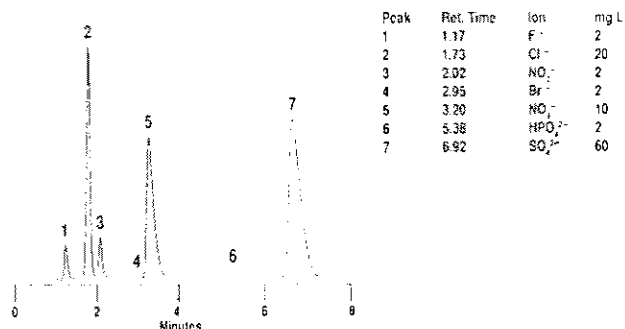


FIG. 1 Chromatogram Showing Separation Using the AS4A Column

TABLE 1 Approximate Single Laboratory Detection Limits in Reagent Water^{A,B}

Analyte	Peak No.	Retention Time, min	MDL mg/L
Fluoride	1	1.2	0.01
Chloride	2	1.7	0.02
Nitrite-N	3	2.0	0.004
Bromide	4	2.9	0.01
Nitrate-N	5	3.2	0.002
o-Phosphate	6	5.4	0.003
Sulfate	7	6.9	0.02

^A Data provided by US EPA/EMSL Laboratory, Cincinnati, OH.

^B Column: as specified in 7.1.4.
 Detector: as specified in 7.1.6.
 Eluent: as specified in 8.3.
 Pump rate: 2.0 mL/min.
 Sample loop: 50 µL.

to separate the ions of interest according to their retention characteristics prior to their detection.

4. Summary of Test Method

4.1 An aliquot of sample is injected into an ion chromatograph. The sample is pumped through two columns and a suppressor device and into a conductivity detector. The analytical column and the guard column are packed with an anion exchange resin. Ions are separated based on their affinity for the exchange sites of the resin. The suppressor device contains a cation exchanger which may be regenerated by a flow of dilute acid. The suppressor device reduces the background conductivity of the eluent to a low or negligible level by replacing the cations with the hydrogen ion, thereby converting the anions in the sample to their corresponding acids. The separated anions in their acid form are measured using an electrical-conductivity cell. Anions are identified based on their retention times compared to known standards. Quantitation is accomplished by measuring the peak height or area and comparing it to a calibration curve generated from known standards.

5. Significance and Use

5.1 Ion chromatography provides for both qualitative and quantitative determination of seven common anions, F⁻, Cl⁻, NO₂⁻, HPO₄²⁻, Br⁻, NO₃⁻, and SO₄²⁻, from a single analytical operation requiring only a few millilitres of sample and taking approximately 10 to 15 min for completion.

NOTE 2—This test method may be used to determine fluoride if its peak is in the water dip by adding one mL of eluent (at 100× the

concentration in 8.3) to all 100-mL volumes of samples and standards to negate the effect of the water dip. (See 6.3, and also see 6.4.) The quantitation of unretained peaks should be avoided. Anions such as low molecular weight organic acids (formate, acetate, propionate, etc.) which are conductive may coelute with fluoride and would bias fluoride quantitation in some drinking waters and most waste waters. Alternate analytical columns may not have this limitation. See Appendix X2.

5.2 Anion combinations such as Cl⁻/Br⁻ and NO₂⁻/NO₃⁻, which may be difficult to distinguish by other analytical methods, are readily separated by ion chromatography.

6. Interferences

6.1 Since chloride and nitrite elute very close together, they are potential interferences for each other. It is advisable not to have one of these anions present in a ten-fold excess over the other; that is, Cl⁻/NO₂⁻ ratios higher than 1:10 or 10:1 if both ions are to be quantitated. Alternate columns may provide better separation and thus will allow higher ratios without interference.

6.2 As with other types of chromatography, if one of the sample components is present at very high levels, it may interfere by causing a very large peak on the chromatogram which could mask other peaks present. This type of interference is normally minimized by dilution of the sample (see Annex A1) and in some instances may be corrected if the concentration of that anion is of interest. However, care should be taken not to dilute the analyte concentration below its detectable limit.

6.3 Water from the sample injection will cause a negative peak or dip in the chromatogram when it elutes, because its conductivity is less than that of the suppressed eluent. This dip usually occurs before Cl⁻. Any peak of interest eluting near the water dip must be sufficiently resolved from the dip to be accurately quantitated. Some suggested techniques for elimination of the water dip are described in Appendix X1.

6.4 Due to the effect of the water dip and the interference of organic acids and due to the presence of carbonate ions in the analytical column, the user of this test method is urged to use caution when determining fluoride (see Note 2). If the user wishes to be certain of good results and has interfering anions present when determining fluoride, the eluent can be diluted until separation of fluoride and the interfering ions are accomplished. This will cause an increase in retention time for all anions.

7. Apparatus

7.1 *Ion Chromatograph*—The ion chromatograph should have the following components assembled, as shown in Fig. 2:⁴

7.1.1 *Eluent and Regenerant Containers.*

7.1.2 *Eluent Pump*, capable of delivering 1 to 3 mL/min of eluent at a pressure of up to 2000 psig.

7.1.3 *Guard Column*—Anion exchange column, typically of the same anion exchange material used in the analytical column. The purpose of this column is to protect the

⁴ Available from Dionex Corp., 1228 Titan Way, Sunnyvale, CA 94086. An equivalent may be used. Other manufacturers' components may provide equivalent data.

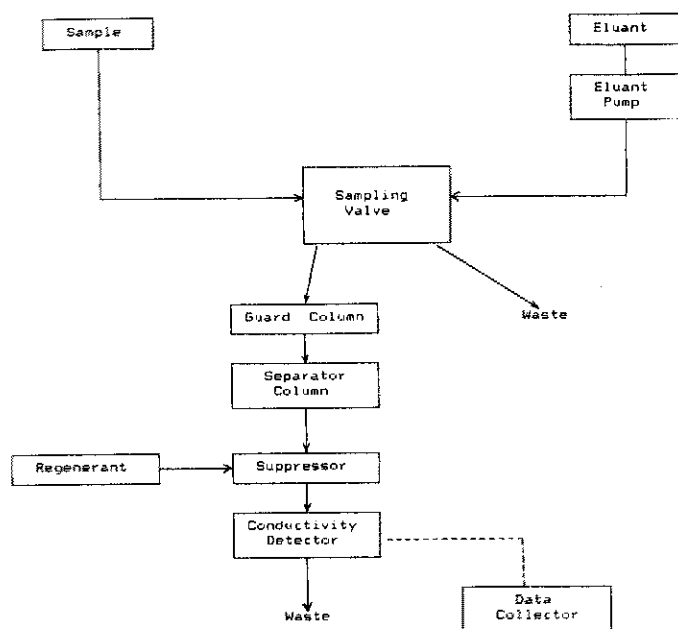


FIG. 2 Schematic of an Ion Chromatograph

analytical column from particulate matter and irreversibly retained materials.

7.1.4 Analytical Column—Anion exchange column capable of separating chloride from the injection void volume, as well as resolving the anions chloride through sulfate.

NOTE 3—Any analytical column may be used. However, the user should be able to achieve the resolution and separation as shown in Fig. 1. See Appendix X2 for an example of an alternate column.

7.1.5 Suppressor Device—A cation exchanger which is regenerated by a flow of dilute acid.

7.1.6 Detector—A low-volume, flow through, temperature-compensated or temperature - controlled electrical conductivity cell equipped with a meter capable of reading from 0 to 1000 $\mu\text{S}/\text{cm}$ on a linear scale.

7.1.7 Recorder, Integrator, Computer—A device compatible with the detector output capable of recording detector response as a function of time for the purpose of measuring peak height or area.

7.1.8 Data System—An electronic integrator, such as is used with gas and liquid chromatographs, may be used to quantitate peak area, as well as peak height. The peak area data can be used in the same way peak height is used to quantitate results. Alternately, a computer with suitable chromatography software may be used.

7.1.9 Sample Loop—A loop on the injection valve that is designed to contain an exact amount of the sample. A common size is 100 μL . The sample volume injected onto the analytical column is controlled by this loop. Use of a larger size loop will usually cause peak broadening and a loop size greater than 1 mL may result in column overloading and nonlinear response. The chromatogram in Fig. 1 uses a 100- μL size sample loop.

7.1.9.1 When injections of volumes larger than the sample loop size are made, any volume above the sample loop size goes to waste. It is considered good technique to flush the sample loop upon injection by injecting 2 to 3 times the sample loop volume.

8. Reagents

8.1 Purity of Reagents—Reagent grade chemicals shall be used in all tests. Unless otherwise indicated, it is intended that all reagents shall conform to the specifications of the Committee on Analytical Reagents of the American Chemical Society, where such specifications are available.⁵ Other grades may be used, provided it is first ascertained that the reagent is of sufficiently high purity to permit its use without lessening the accuracy of the determination.

8.2 Purity of Water—Unless otherwise indicated, references to water shall be understood to mean reagent water conforming to Specification D 1193, Type II. Column life may be extended by passing Type II water through a 0.22- μm filter prior to use. Freshly prepared water should be used for making the standards intended for calibration. The detection limits of this test method will be limited by the purity of the water and reagents used to make the standards. The purity of the water may be checked by use of this test method. Anion concentrations of less than 0.2 $\mu\text{g}/\text{L}$ each are typical of this type of water.

8.3 Eluent—Dissolve 0.2856 g of sodium bicarbonate (1.7 mM) and 0.3816 g of sodium carbonate (1.8 mM) in water and dilute to 2 L with water. The eluent should be degassed prior to use. Other eluents may also prove to be acceptable, provided they give the proper resolution between the component peaks. This eluent will act as a growth media for algae. For this reason the eluent should not be kept for longer than one month.

NOTE 4—Use of other eluents may change the order of elution of the anions from that using the carbonate-bicarbonate eluent.

8.4 Suppressor Regenerant Solution—Cautiously add 3 mL of H_2SO_4 (sp gr 1.84) to 4 L of water. Suppressors using electrochemical generation of acid will not require this regenerant.

8.5 Stock Solutions:

8.5.1 Bromide Stock Solution (1.00 mL = 1.00 mg Br^-)—Dry approximately 2 g of sodium bromide (NaBr) for 6 h at 150°C and cool in a desiccator. Dissolve 1.2877 g of the dried salt in water and dilute to 1 L with water.

8.5.2 Chloride Stock Solution (1.00 mL = 1.00 mg Cl^-)—Dry sodium chloride (NaCl) for 1 h at 100°C and cool in a desiccator. Dissolve 1.648 g of the dry salt in water and dilute to 1 L with water.

8.5.3 Fluoride Stock Solution (1.00 mL = 1.00 mg F^-)—Dissolve 2.210 g of sodium fluoride (NaF) in water and dilute to 1 L with water.

8.5.4 Nitrate Stock Solution (1.00 mL = 1.00 mg NO_3^-)—Dry approximately 2 g of sodium nitrate (NaNO_3) at 105°C for 48 h. Dissolve exactly 1.371 g of the dried salt in water and dilute to 1 L with water.

8.5.5 Nitrite Stock Solution (1.00 mL = 1.00 mg NO_2^-)—Place approximately 2 g of sodium nitrite (NaNO_2) in a 125-mL beaker and dry to constant weight (about 24 h) in a desiccator containing concentrated H_2SO_4 . Dissolve 1.500 g of the dried salt in water and dilute to 1 L with water. Store

⁵ "Reagent Chemicals, American Chemical Society Specifications," Am. Chem. Soc., Washington, DC. For suggestions on the testing of reagents not listed by the American Chemical Society, see "Analytical Standards for Laboratory Chemicals," BDH Ltd., Poole, Dorset, U.K., and the "United States Pharmacopeia."

in a sterilized glass bottle. Refrigerate and prepare monthly.

NOTE 5—Nitrite is easily oxidized, especially in the presence of moisture, and only fresh reagents are to be used.

NOTE 6—Prepare sterile bottles for storing nitrite solutions by heating for 1 h at 170°C in an air oven.

8.5.6 Phosphate Stock Solution (1.00 mL = 1.00 mg HPO_4^{-2})—Dissolve 1.433 g of potassium dihydrogen phosphate (KH_2PO_4) in water and dilute to 1 L with water.

8.5.7 Sulfate Stock Solution (1.00 mL = 1.00 mg SO_4^{-2})—Dry sodium sulfate (Na_2SO_4) for 1 h at 105°C and cool in a desiccator. Dissolve 1.479 g of the dried salt in water and dilute to 1 L with water.

8.6 Anion Working Solutions—Prepare a blank and at least 3 different working standards containing the anions of interest. The combination anion solutions should be prepared in volumetric flasks. These standards must be prepared fresh daily. The concentration range for the three standards will be dependent on the levels expected in the samples. If desired, a single standard may be prepared that contains all six anions.

8.6.1 The user should select the ranges of the three standards so as to cover the entire range of sample concentration expected. The ranges chosen should all fall into one attenuation setting. If a second attenuation setting must be used, it must be calibrated using three standards and a blank. The standard concentrations given in Tables 2 and 3 are for example purposes.

9. Sampling

9.1 Collect the sample in accordance with Practices D 1066 and D 3370 as applicable.

9.2 Analyze the samples as soon as possible after collection. Preservation by refrigeration at 4°C is required for nitrite, nitrate, or phosphate.

9.3 Filter the samples containing particulates through a prewashed 0.22- μm filter prior to analysis to avoid fouling or clogging the resin of the columns.

10. Calibration

10.1 Determination of Retention Times:

10.1.1 The retention time for each anion is determined by injecting a standard solution containing only the anion of interest and noting the time required for a peak to appear on the chromatogram. Retention times vary with operating conditions and are influenced by the concentration of ion(s)

TABLE 2 Preparation of Standard Solutions for Instrument Calibration

Anion	High-Range Standard		Intermediate-Range Standard, mg/L	Low-Range Standard, mg/L
	Millilitres of Each Stock Solution (1.00 mL = 1.00 mg), Diluted to 1000 mL	Anion Concentration, mg/L		
Fluoride (F^-)	10	10	1.0	0.2
Chloride (Cl^-)	10	10	1.0	0.2
Nitrite (NO_2^-)	20	20	2.0	0.4
Phosphate (HPO_4^{-2})	50	50	5.0	1.0
Bromide (Br^-)	10	10	1.0	0.2
Nitrate (NO_3^-)	30	30	3.0	0.6
Sulfate (SO_4^{-2})	100	100	10.0	2.0

TABLE 3 Preparation of Standard Solutions for Determination of Retention Times

Stock Solution (1 mL = 1.00 mg)	Volume of Stock Solution per Litre of Water, mL	Anion Concentration, mg/L
Fluoride	4	4
Chloride	4	4
Nitrite	10	10
Phosphate	50	50
Bromide	10	10
Nitrate	30	30
Sulfate	50	50

present. Prepare separate standard solutions in accordance with Table 3 by pipetting the designated aliquots of stock solutions prepared in Section 8 (8.5.1 through 8.5.7) into separate 1-L volumetric flasks. Analyze each standard of interest as defined in Section 11. Note the time in minutes for each peak to appear on the chromatogram.

NOTE 7—Some operators have reported unusually large shifts in retention time for nitrate with changes in concentration. If this occurs, care must be taken to ensure integration of the correct peak when integration is used for calculation.

10.1.2 Concentrations other than those listed in Table 3 may be used if they better approximate concentrations expected in the samples.

10.1.3 Retention times as well as elution order vary with the column used. See Fig. 1 for example elution orders.

10.2 Analyze the blank and each of the prepared calibration solutions described in 8.7 in accordance with the defined procedure (see Section 11).

NOTE 8—If the concentrations of the sample ions of interest are known or estimated, the concentration of standard solutions prepared for instrument calibration may be varied to better approximate or bracket the concentration range of interest. Anions of no interest may be omitted.

NOTE 9—The mid-range combination anion standard may be used to verify resolution of all seven anions.

NOTE 10—Each analytical curve should be established using only one scale setting. Changing the scale setting may result in a slight change in the slope of the analytical curve.

10.3 Prepare analytical curves for each anion of interest by plotting on linear graph paper peak height or peak area versus the nominal concentration of the anion standard solution. Electronic integrators and computer based chromatography data systems provide this function automatically.

NOTE 11—Some operators have reported a shift in slope of the phosphate calibration curve at approximately 30 mg/L HPO_4^{-2} . If such a shift in slope is observed, additional standard solutions covering the entire range of concentration should be prepared and analyzed in order to accurately define the slope of the curve. If an integrator is being used, it may be necessary to manually calculate phosphate concentrations above 30 mg/L HPO_4^{-2} in order to obtain maximum accuracy.

11. Procedure

11.1 Set up the ion chromatograph according to the manufacturer's instructions.

11.1.1 The detector ranges are variable. Normal operating ranges are from 1 to 30 $\mu\text{S}/\text{cm}$ full scale. The range setting required for analyses will depend on the concentration of

TABLE 4 Determination of Bias for Fluoride

Water	Amount Added, mg/L	Amount Found, mg/L	S_i	S_o	Bias, %
Reagent	0.26	0.25	0.08	0.11	-3.8
	0.34	0.29	0.11		-14.7
	2.12	2.12	0.07	0.12	0.0
	2.55	2.48	0.14		-2.7
	6.79	6.76	0.20	0.19	-0.4
	8.49	8.46	0.30		-0.4
Drinking	0.26	0.24	0.08	0.05	-7.7
	0.34	0.34	0.11		0.0
	2.12	2.09	0.18	0.06	-1.4
	2.55	2.55	0.16		0.0
	6.79	6.84	0.54	0.25	+0.7
	8.49	8.37	0.75		-1.4
Waste	0.26	0.25	0.15	0.06	-3.8
	0.34	0.32	0.08		-5.9
	2.12	2.13	0.22	0.15	+0.5
	2.55	2.48	0.16		-2.7
	6.79	6.65	0.41	0.20	-2.1
	8.49	8.27	0.36		-2.6

TABLE 5 Determination of Bias for Chloride

Water	Amount Added, mg/L	Amount Found, mg/L	S_i	S_o	Bias, %
Reagent	0.78	0.79	0.17	0.29	+1.3
	1.04	1.12	0.46		+7.7
	6.50	6.31	0.27	0.14	-2.9
	7.80	7.76	0.39		-0.5
	20.8	20.7	0.54	0.62	-0.5
	26.0	25.9	0.58		-0.4
Drinking	0.78	0.54	0.35	0.20	-30.8
	1.04	0.51	0.38		-51.0
	6.50	5.24	1.35	1.48	-19.4
	7.80	6.02	1.90		-22.8
	20.8	20.0	2.26	1.14	-3.8
	26.0	24.0	2.65		-7.7
Waste	0.78	0.43	0.32	0.39	-44.9
	1.04	0.65	0.48		-37.5
	6.50	4.59	1.82	0.83	-29.4
	7.80	5.45	2.02		-30.1
	20.8	18.3	2.41	1.57	-11.8
	26.0	23.0	2.50		-11.5

TABLE 6 Determination of Bias for Nitrite-Nitrogen

Water	Amount Added, mg/L	Amount Found, mg/L	S_i	S_o	Bias, %
Reagent	0.36	0.37	0.04	0.04	+2.8
	0.48	0.48	0.06		0.0
	3.00	3.18	0.12	0.06	+6.0
	3.60	3.83	0.12		+6.4
	9.60	9.84	0.36	0.26	+2.5
	12.0	12.1	0.27		+0.6
Drinking	0.36	0.30	0.13	0.03	-16.7
	0.48	0.40	0.14		-16.7
	3.00	3.02	0.23	0.12	+0.7
	3.60	3.62	0.22		+0.6
	9.60	9.59	0.44	0.28	-0.1
	12.0	11.6	0.59		-3.1
Waste	0.36	0.34	0.06	0.04	-5.6
	0.48	0.46	0.07		-4.2
	3.00	3.18	0.13	0.10	+6.0
	3.60	3.76	0.18		+4.4
	9.60	9.74	0.49	0.26	+1.5
	12.0	12.0	0.56		+0.3

TABLE 7 Determination of Bias for Bromide

Water	Amount Added, mg/L	Amount Found, mg/L	S_i	S_o	Bias, %
Reagent	0.63	0.69	0.11	0.05	+9.5
	0.84	0.85	0.12		+1.2
	5.24	5.21	0.22	0.21	-0.6
	6.29	6.17	0.35		-1.9
	16.8	17.1	0.70	0.36	+1.6
	21.0	21.3	0.93		+1.5
Drinking	0.63	0.63	0.13	0.04	0.0
	0.84	0.81	0.13		-3.6
	5.24	5.11	0.23	0.13	-2.5
	6.29	6.18	0.30		-1.7
	16.8	17.0	0.55	0.57	+0.9
	21.0	20.9	0.65		-0.4
Waste	0.63	0.63	0.15	0.09	0.0
	0.84	0.85	0.15		+1.2
	5.24	5.23	0.36	0.11	-0.2
	6.29	6.27	0.46		-0.3
	16.8	16.6	0.69	0.43	-1.0
	21.0	21.1	0.63		+0.3

anions in the sample and should be chosen accordingly.

11.2 Equilibrate the system by pumping eluent through the analytical column and suppressor device and the detector until a stable baseline is obtained (approximately 15 to 20 min) and approximately 10 to 15 $\mu\text{S}/\text{cm}$ background conductivity for the recommended eluent. This equilibration normally can be accomplished while the samples and standards are being prepared.

11.3 Load 2 to 3 mL of sample into the sample entry port using a syringe. Inject the sample into the eluent stream and record the ion chromatogram (see Fig. 1).

NOTE 12—Most plastic disposable syringes, if used, come in packages labeled as "sterile," but this does not necessarily imply that they are free of anions; therefore, these syringes should be flushed with water and sample or standard before use to minimize contamination.

12. Calculation

12.1 Compare the peak heights or areas noted for the anion(s) in the samples to the calibration curves prepared in 10.3 to determine the anion concentration in milligrams per litre:

$$\text{Anion concentration, mg/L} = A \times F$$

where:

A = milligrams per litre read from appropriate calibration, and

F = dilution factor if sample was diluted prior to analyses.

13. Precision and Bias⁶

13.1 The precision and bias data presented in this test method meet the requirements of the 1986 edition of Practice D 2777.

13.2 The collaborative test of this test method was performed in reagent water, drinking water, and a waste water of choice by 19 laboratories using one operator each. Six levels of concentration were used for seven anions, producing three youden pairs. Each youden pair was used to calculate the single operator precision (S_o).

13.3 The precision and bias of this test method for each

⁶ Supporting data are available from ASTM Headquarters. Request RR: D19-1147.

TABLE 8 Determination of Bias for Nitrate-Nitrogen

Water	Amount Added, mg/L	Amount Found, mg/L	S_t	S_o	Bias, %
Reagent	0.42	0.42	0.04	0.02	0.0
	0.56	0.56	0.06		0.0
	3.51	3.34	0.15	0.08	-4.8
	4.21	4.05	0.28		-3.8
	11.2	11.1	0.47	0.34	-1.1
	14.0	14.4	0.61		+2.6
Drinking	0.42	0.46	0.08	0.03	+9.5
	0.56	0.58	0.09		+3.6
	3.51	3.45	0.27	0.10	-1.7
	4.21	4.21	0.38		0.0
	11.2	11.5	0.50	0.48	+2.3
	14.0	14.2	0.70		+1.6
Waste	0.42	0.36	0.07	0.06	-14.6
	0.56	0.40	0.16		-28.6
	3.51	3.19	0.31	0.07	-9.1
	4.21	3.84	0.28		-8.8
	11.2	10.9	0.35	0.51	-3.0
	14.0	14.1	0.74		+0.4

TABLE 9 Determination of Bias for *ortho*-Phosphate

Water	Amount Added, mg/L	Amount Found, mg/L	S_t	S_o	Bias, %
Reagent	0.69	0.69	0.06	0.06	0.0
	0.92	0.98	0.15		+6.5
	5.77	5.72	0.36	0.18	-0.9
	6.92	6.78	0.42		-2.0
	18.4	18.8	1.04	0.63	+2.1
	23.1	23.2	0.35		+0.4
Drinking	0.69	0.70	0.17	0.17	+1.4
	0.92	0.96	0.20		+4.3
	5.77	5.43	0.52	0.40	-5.9
	6.92	6.29	0.72		-9.1
	18.4	18.0	0.68	0.59	-2.2
	23.1	22.6	1.07		-2.0
Waste	0.68	0.64	0.26	0.09	-7.2
	0.92	0.82	0.28		-10.9
	5.77	5.18	0.66	0.34	-10.2
	6.92	6.24	0.74		-9.8
	18.4	17.6	2.08	1.27	-4.1
	23.1	22.4	0.87		-3.0

anion for reagent water, drinking water, and waste water are shown in Table 4 through 10.

13.4 Some of the bias statements, for example chlorine and sulfate, may be misleading due to spiking of small increments of the anion onto large naturally occurring concentrations of the same anion.

13.5 All data in Tables 4 through 10 were obtained using a non-metallic pump surface so as to minimize potential metallic contamination to the analytical columns.

TABLE 10 Determination of Bias for Sulfate

Water	Amount Added, mg/L	Amount Found, mg/L	S_t	S_o	Bias, %
Reagent	2.85	2.83	0.32	0.52	-0.7
	3.80	3.83	0.92		+0.8
	23.8	24.0	1.67	0.68	+0.8
	28.5	28.5	1.56		-0.1
	76.0	76.8	3.42	2.33	+1.1
	95.0	95.7	3.59		+0.7
Drinking	2.85	1.12	0.37	0.41	-60.7
	3.80	2.26	0.97		-40.3
	23.8	21.8	1.26	0.51	-8.4
	28.5	25.9	2.48		-9.1
	76.0	74.5	4.63	2.70	-2.0
	95.0	92.3	5.19		-2.8
Waste	2.85	1.89	0.37	0.24	-33.7
	3.80	2.10	1.25		-44.7
	23.8	20.3	3.19	0.58	-14.7
	28.5	24.5	3.24		-14.0
	76.0	71.4	5.65	3.39	-6.1
	95.0	90.3	6.80		-5.0

TABLE 11 Quality Control (QC) Sample Precision and Bias Data from Linear Regression Equations in mg/L

Water	Anion	True Value of QC's	Mean Rec. of QC's	S_t of QC's
Reagent	F ⁻	2.01	1.98	0.128
	Cl ⁻	10.00	9.94	0.490
	NO ₂ ⁻	5.00	5.18	0.170
	Br ⁻	5.00	5.02	0.265
	NO ₃ ⁻	5.01	4.93	0.256
	HPO ₄ ⁻	7.01	7.03	0.383
	SO ₄ ⁻²	25.00	25.17	1.58
Drinking	F ⁻	2.01	2.00	0.203
	Cl ⁻	10.00	8.55	1.88
	NO ₂ ⁻	5.00	4.97	0.295
	Br ⁻	5.00	4.95	0.240
	NO ₃ ⁻	5.01	5.03	0.336
	HPO ₄ ⁻	7.01	6.71	0.540
	SO ₄ ⁻²	25.00	23.07	2.02
Waste	F ⁻	2.01	1.97	0.170
	Cl ⁻	10.00	7.95	2.33
	NO ₂ ⁻	5.00	5.15	2.35
	Br ⁻	5.00	5.00	0.320
	NO ₃ ⁻	5.01	4.72	0.356
	HPO ₄ ⁻	7.01	6.53	0.746
	SO ₄ ⁻²	25.00	21.98	3.13

13.6 Table 2 compares the mean recoveries for each ion obtained by all laboratories with their true values. This table also shows the overall (S_t) precision that was achieved by the participating laboratories.

14. Keywords

14.1 anions; drinking water; ion chromatography; reagent water; waste water

ANNEX

(Mandatory Information)

A1. DETERMINATION OF RANGE OF TEST METHOD

A1.1 The upper range of the test method for an injected aliquot varies for each sample type and is limited by the capacity of the analytical column. If the ionic strength of the sample to be analyzed exceeds column capacity, the ions of interest will not be separated as shown in Fig. 1. Retention times will also decrease as the ionic strength of the sample increases. The upper range of the test method for the ion(s) of interest in the sample matrix may be determined by the following procedure. The procedure is defined for currently available analytical columns supplied by the appropriate manufacturer.

A1.1.1 Bypass the analytical column (only the suppressor column is on-line) and inject 100 μ L of 5.8 g/L (0.01 meq)

solution of sodium chloride (NaCl). This represents 10 % of the analytical column capacity, that is, 0.1 meq. Note the conductivity response.

A1.1.2 Inject 100 μ L of the sample of interest (analytical column bypassed). Note the conductivity response.

A1.1.3 If the conductivity response noted in A1.1.2 is less than or equal to the conductivity response noted in the sodium chloride solution in A1.1.1, column capacity has not been exceeded.

A1.1.4 If the conductivity response noted in A1.1.2 is greater than the conductivity response noted for the sodium chloride solution in A1.1.1, the test sample must be diluted prior to injection into the ion chromatograph.

APPENDIX

(Nonmandatory Information)

X1. TECHNIQUES FOR ELIMINATING WATER DIP

X1.1 A method of eliminating the water dip is to introduce concentrations of carbonate and bicarbonate into the sample that closely approximate that of the eluent used for analysis. Adjustment of sample background may be accomplished in two different ways:

X1.1.1 Dilute the sample with eluent if sample dilution is required prior to analysis.

X1.1.2 Add an equivalent of 1.0 mL of a prepared eluent concentrate (solution that is 100 times more concentrated than the eluent used for analysis) per 100 mL of sample.

X1.2 Standard solutions must be prepared as directed in X1.1.2. It is important to prepare a blank using reagent water and eluent (100:1) to compensate for any anionic impurities present.

X2. ALTERNATE ANALYTICAL COLUMN

X2.1 In time, new ion exchange analytical columns become available which overcome some of the limitations of older columns. Figure X2.1 shows an example of an analytical column⁷ which has the following improvements over the column shown in Fig. 1. This column retains fluoride sufficiently to separate it from the void volume (water dip), while also separating it from low molecular organic acids, such as acetate. In addition, it gives better separation of nitrite from chloride. This column also provides these improvements, while still eluting sulfate in under 14 min.

⁷ Available from Dionex Corp., 1228 Titan Way, Sunnyvale, CA 94086. Other manufacturer's columns may provide equivalent performance.

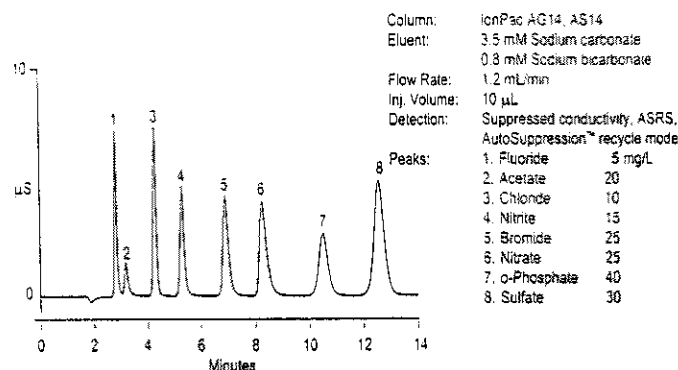


FIG. X2.1 Chromatogram Showing Separation Using the AS14 Column

The American Society for Testing and Materials takes no position respecting the validity of any patent rights asserted in connection with any item mentioned in this standard. Users of this standard are expressly advised that determination of the validity of any such patent rights, and the risk of infringement of such rights, are entirely their own responsibility.

This standard is subject to revision at any time by the responsible technical committee and must be reviewed every five years and if not revised, either reapproved or withdrawn. Your comments are invited either for revision of this standard or for additional standards and should be addressed to ASTM Headquarters. Your comments will receive careful consideration at a meeting of the responsible technical committee, which you may attend. If you feel that your comments have not received a fair hearing you should make your views known to the ASTM Committee on Standards, 100 Barr Harbor Drive, West Conshohocken, PA 19428.